

CERVMAX™Tx

Cervical Orthosis Collar
with Thoracic Extension

Benefits

- » New thoracic extension for maximum stability
- » Uncompromising comfort and support
- » CoolMax® liner
Cooler than any other fabric
- » New reimbursement code
- » Updated sizes



Model 0000



Model 0000



Model 0000




Model 0000

1605 Cervmax™Tx PRODUCT INFORMATION

CERVMAX™Tx
Cervical Orthosis Collar
with Thoracic Extension

PMT offers the new Cervmax™Tx for physicians who desire a cervical collar with an extension for additional thoracic support.

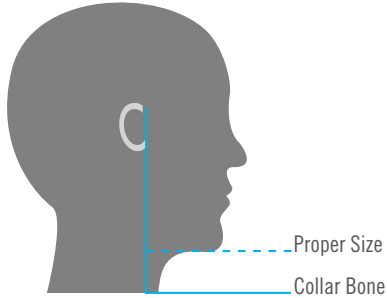

SIZE	SHOULDER TO CHIN
Short	1 1/2" - 2 1/4"
Stout	2 1/4" - 2 3/4"
Small	2 1/4" - 2 3/4"
Medium	2 3/4" - 3 1/4"
Large	3 1/4" - 3 3/4"
Tall	3 3/4" - 4 1/4"

The unique V-Back design of the Cervmax™Tx Collar provides maximum occipital support, cradling the back of the head and is flexible enough to conform to most neck shapes. Cervmax™Tx features less shell coverage in the chin area which reduces skin irritation and maintains the maximum mandibular support to control tilt (coronal flexion and extension). The slotted front portion aids in even weight distribution across the sternal notch area, increasing patient comfort with less pressure. Lightweight materials are MRI compatible, silicone-free and latex-free.

Supplied With Cervmax™Tx

Each Cervmax™Tx Collar is provided standard with an extra set of replacement pads and size identification card.

Our Commitment to Customer Satisfaction

PMT is committed to providing you the highest quality products and the finest customer service. Local sales representatives are available to serve you.



Cervmax™Tx Features	Cervmax™Tx Benefits
A. Patented V-Back Design	Provides superior occipital support and immobilization.
B. Thoracic Extension	Provides additional thoracic support to the patient.
C. CoolMax® Liner	Wicks moisture away from skin, allows skin to breathe. The hook and loop system allows for easy removal of washable, fast-drying, durable liner.
D. Velcro™ Straps	Provides additional thoracic support to the patient.
E. Slotted Shell	Allows airflow and aids in breathability.

▶ **Local sales representatives are available 24 hours a day, 7 days a week, 365 days a year to assist you. Contact PMT® Corporation for more information.**

These products are protected by one or more of the following U.S. patents: 6,872,188 & 6,494,854.

PMT and Mankind are registered trademarks and Cervmax™Tx is a trademark of PMT Corporation. CoolMax® is a DuPont trademark. © 2010 PMT Corporation. All Rights Reserved

FDA REGISTERED AND APPROVED 

PMT CORPORATION
QUALITY, INNOVATION AND CUSTOMER DRIVEN

PMT CORPORATION
1500 Park Road
Chanhassen, MN USA 55317
www.pmtcorp.com

Call Toll Free ... 24 Hours a Day/365 Day a Year
1-800-MANKIND® (626-5463)
Phone 952-470-0866 » Fax 952-470-0865